

Pergoveris™

(follitropin alfa and lutropin alfa for injection)



Prefilled Pergoveris Pen

How to use your Pergoveris **pre-filled** pen

Only 3 steps



1

Attaching the needle

- Remove the pen cap.
- Make sure that the dose feedback window is set at '0'.
- Remove the peel-off seal from the outer needle cap.
- Screw the threaded tip of the pen into the outer needle cap until you feel slight resistance. Take care not to attach it too tightly as it may become difficult to remove the needle after the injection.



2

Setting the dose

- Turn the dose-setting knob until the dose prescribed by your doctor shows in the dose feedback window. You can adjust the dose in 12.5 IU increments.
- If needed, you can correct the dose by turning the knob back anti-clockwise.
- Always check that the dose feedback window shows the correct dose before injecting.



3

Injecting your dose

- Slowly push the needle into the skin at a 90° angle.
- Push down on the dose-setting knob as far as it will go. The larger the dose, the longer the injection will take. Hold the dose-setting knob down for at least 5 seconds. The dose feedback window will return to '0'.
- Pull the needle out of the skin, while keeping the dose-setting knob pressed down.
- Choosing a different injection site each day will help to minimize skin irritation.

Pergoveris™

(follitropin alfa and lutropin alfa for injection)




Prefilled Pergoveris Pen

Safety

Always dispose of needles safely as recommended by your healthcare professional. If needed, ask your doctor or pharmacist for more information.

NEVER REUSE NEEDLES

Storage:

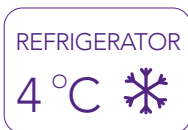

Please write the date of the first use of the pen.
Once opened, the pre-filled pen may be stored for a maximum of 28 days outside of the refrigerator at room temperature (at 25°C).

Do not use the product after the expiry date (EXP).
The expiry date refers to the last day of that month.
Do not use the product if you notice that the liquid contains particles or is not clear.

Merck Oy
Keilaranta 6, 02150 Espoo
www.merck.fi

FI/PER-00003

MERCK